



**INVACARE 20" X 18"H BACKREST (ANGLED BENDS)**



**WELD NUTS AND SCREWS**

*INVACARE 20" X 18"H BACKREST WITH BODYPOINT FRONT-PULL TOP STRAPS AND LOWER STRAPS FITTED USING STANDARD SCREWS AND NUTS (1/4" OR M6). WELD NUTS MAY BE SUBSTITUTED FOR A LOWER-PROFILE ON THE INNER SURFACE TO AVOID DAMAGE TO UPHOLSTERY.*


|      |   |                                     |                            |                  |
|------|---|-------------------------------------|----------------------------|------------------|
| note | description                               |                                     |                            |                  |
|      | <b>MODIFICATION TO INVACARE BACKRESTS</b> |                                     |                            |                  |
|      | drawing no.<br><i>D11-0111-01</i>         | date originated<br><i>1/11/2011</i> | drawn by<br><i>M. KOSH</i> | rev.<br><b>A</b> |



**Bodypoint**<sup>®</sup>

technical bulletin

FOR FRONT-PULL HARNESS TOP STRAPS

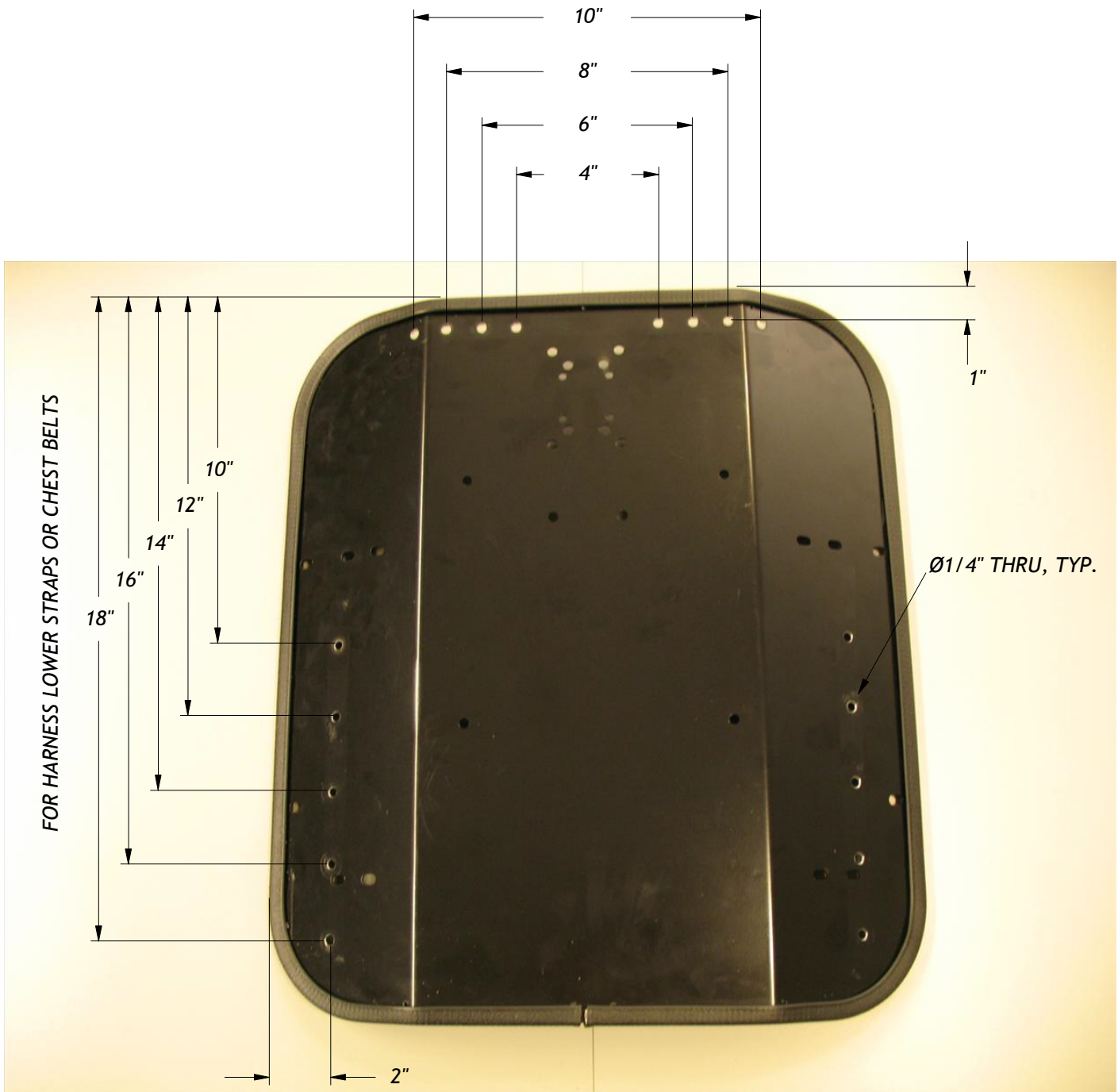

**INVACARE 20" X 18"H BACKREST (ANGLED BENDS)  
CURVED BACKRESTS OF THE SAME SIZE CAN USE  
THIS SAME HOLE LAYOUT**

note

description

**MODIFICATION TO INVACARE BACKRESTS**

drawing no.

D11-0111-01

date originated

1/11/2011

drawn by

M. KOSH



rev.

A



LOWER STRAPS CAN BE FITTED TO EXISTING SLOTS USING WELD NUTS

FOR FRONT-PULL HARNESS TOP STRAPS

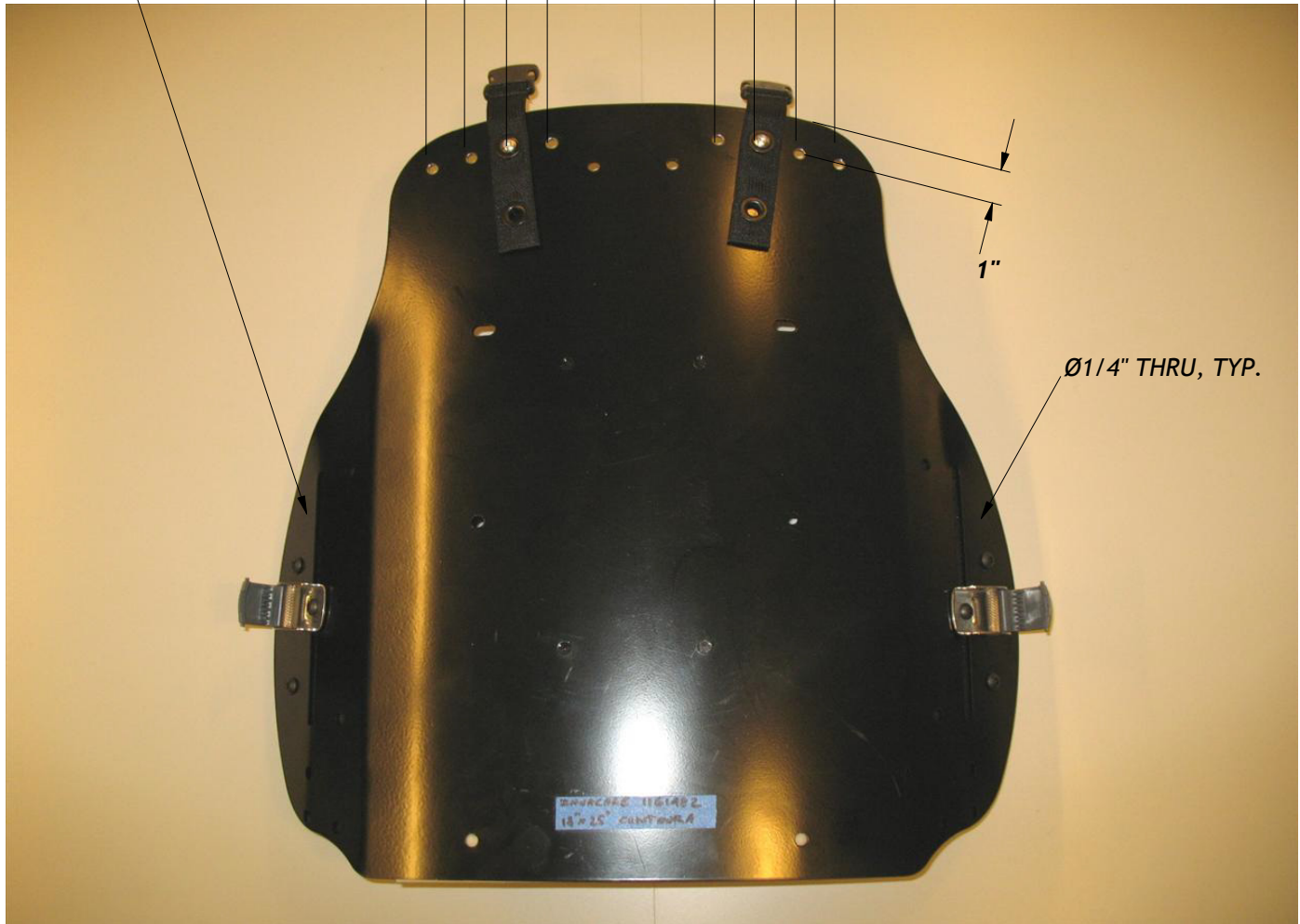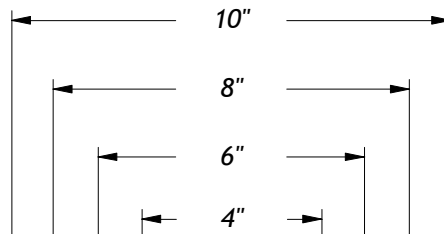

Ø1/4" THRU, TYP.

INVACARE CONTOUR BACKREST HAS EXISTING SLOTS ALONG EACH SIDE WHICH CAN BE USED AS MOUNTING LOCATIONS FOR THE LOWER STRAPS (IF NO HARDWARE OBSTRUCTIONS EXIST). IF HIGHER OR LOWER MOUNTING LOCATIONS ARE NEEDED FOR CHEST BELTS OR LOWER STRAPS ON LARGER USERS, EXISTING HOLES OR ADDING HOLES MAY BE REQUIRED.

INVACARE 20" X 18"H BACKREST (ANGLED BENDS) CURVED BACKRESTS OF THE SAME SIZE CAN USE THIS SAME HOLE LAYOUT

|      |   |                              |                     |           |
|------|---|------------------------------|---------------------|-----------|
| note | description                               |                              |                     |           |
|      | <b>MODIFICATION TO INVACARE BACKRESTS</b> |                              |                     |           |
|      | drawing no.<br>D11-0111-01                | date originated<br>1/11/2011 | drawn by<br>M. KOSH | rev.<br>A |





**Bodypoint®**  
technical bulletin



*MONOFLEX CHEST BELTS CAN BE FITTED TO EITHER ANY STYLE OF BACKREST USING SCREWS TO ATTACH THE INCLUDED SLOTTED MOUNTING PLATE SHOWN HERE.*



*AT RIGHT, THE INVACARE InTouch PCS BACKREST WITH A MONOFLEX INSTALLED*

note

description

**MODIFICATION TO INVACARE BACKRESTS**

drawing no.

*D11-0111-01*

date originated

*1/11/2011*

drawn by

*M. KOSH*

rev.

*A*